

Installing MagAO-X on the Telescope

This procedure describes how to install the MagAO-X instrument on the Magellan Clay Telescope.

Estimated Time to Complete: 5 hours

This document can be downloaded as a PDF: [Installing MagAO-X on the Telescope](#)

Initial Conditions

- Instrument in LCO cleanroom on legs

Preparations

- Shutdown MagAO-X
- Prepare Legs for Move
 - Turn off table legs (table air power off, PEPS II rocker switch off) and remove air connection from wall
 - Lie under table and remove all three (#1,2,& 3) geophones from bottom of table
 - Place coiled up geophones and cables into basket
 - Unplug PEPS II power from side of rack
 - Remove the taped down exhaust line from the clean room floor
- Remove all cables
 - See detailed procedure for removing 2K DM cables
- Electronics Rack
 - Ensure that roll-out shelves are restrained
 - If not installed, install side panels
 - Close and lock doors
 - Tape keys down
- Instrument
 - Remove eyepiece
 - Place freezer baggie over top periscope mirror so it does not fall on pyramid

- Turn off blower, and remove hose
- Tape over any exposed holes (from cables, etc)
- Secure any loose cables
- Shrink wrap the instrument
- Install solar blanked over shrink wrap
- Cart and Rigging
 - Verify all cart hardware is in-hand
 - Verify two wire harnesses are in hand
 - Partially assemble cart, leaving one long side off
 - Store cart out of the way so table can roll into unpacking area

Rig Onto Cart

- Lower table legs onto the casters by turning the 16 leveling bolts, and remove the 12 metal pads
- Roll the instrument out of the cleanroom


NOTE: At least 4 people required, 2 to push, 2 to hold the cleanroom sides open

CAUTION: only push on the legs, not the table itself

- Continue rolling the instrument through the garage door into the unpacking area
- Move the cart, currently with 3 sides, around the legs

CAUTION: be careful to not bump the legs with the cart

- Attach the 4th side of the cart
- Ensure that the 8 large bolts on the cart are snug but not tight
- Install the guide pins along one side of table

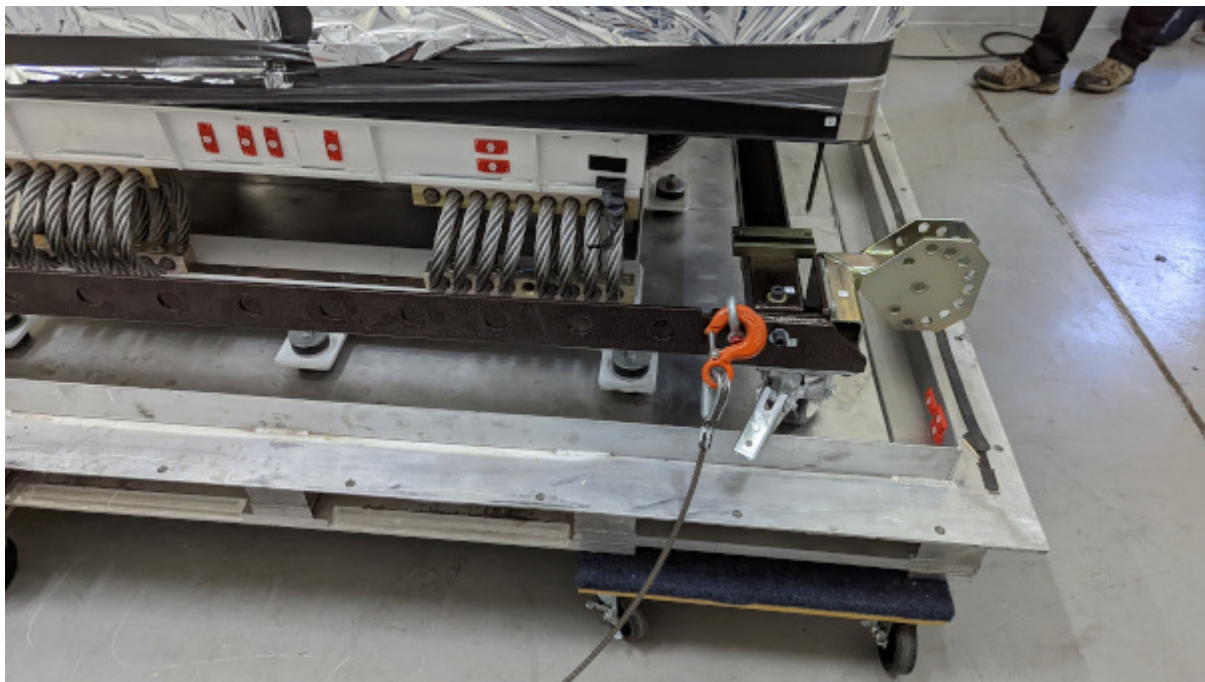

One of the guide pins installed, visible at the lower right corner of table

- Attach the lifting wire-harness to each side of the cart
- Attach the load spreader with straight extensions to the crane, using a crane scale



The load spreader attached to the crane for lifting the cart

- Place the load spreader in the center position (the cart is symmetric)

- Lift the load spreader, and position it over the instrument
- Being careful not to bump the instrument, lower the load spreader and attach the lifting harness D rings. Use 2x shackles to extend the length to reach the cart on the floor.



Lifting harnesses attached with shackle extensions

- Position a person at each corner of the cart

CAUTION: Do not allow the cart to bump the legs or the table uncontrolled

- Slowly lift the cart (**320 lbs**) until it is touching the bottom of the table. Ensure the guide pins go into the outboard slots in the brackets.



The cart being lifted to the bottom of the table.

- Install the 2 bolts opposite the guide-pins. Loosen bolts on the cart as needed.
- Remove the 2 guide-pins, and replace with bolts. There are a total of 4 bolts now installed.
- Once the cart is bolted to the table bottom, while **320 lbs** is still on the crane, tighten all cart bolts. Do not over-tighten, make 1/4 turn after the washers are no longer free. This is to avoid excessive stress on the table.
- Install the triangle stabilizing ropes between the crane hook and the lifting fixture IAW the below figure. Tighten, but do not cause them to pick the load.



The triangle stabilizing ropes should be tight, but not become the lifting point for the load.

CAUTION: be sure that the load spread does not hit the cart when being repositioned.

- Reposition the load spreader center to the instrument + cart position marked on it.
- Ensure that there is room to move the legs out from under the table, opening the garage door into the cleanroom if necessary.
- Position a person at each corner of the cart to stabilize it during the lift.
- Position two people to remove the legs from under the table
- Begin SLOWLY lifting the table off the legs. The weight of the table + cart is **1920 lbs**
- Go up 0.5inch then stop and inspect the platters. Use a long flathead screwdriver to very gently pry any stuck platters off bottom of table.
- Move the legs out from under the table.



The cart and instrument ready to be set down on the wheels, with legs out of the way.

- Set the cart down on its wheels.

Transport MagAO-X To Clay

- Move MagAO-X onto the lift gate (using plates over the gap)



MagAO-X moved onto the pentalift. Note the metal plates across the gap.

- Raise the lift-gate to the height of the flatbed truck



Pentalift raised to truck height.

- Install the come-along and place the plates across the gap.
- Move MagAO-X onto the truck using the come-along
- Secure the instrument by strapping the cart down at 4 points as illustrated in the below figure.



MagAO-X will be strapped to the Isuzu flatbed.

- Slowly drive the truck to the summit.
- Ensure that the lift gate at the summit has been adjusted for slow smooth operation as is done for the asm
- Back the flatbed truck up to the lift gate.



MagAO-X will be unloaded at the telescope using the lift gate, adjusted for slow operation as it is for the ASM.

- Next, using the come-along, carefully move MagAO-X onto the lift gate.



Use the come-along to move MagAO-X on the gate

- Very slowly so as to minimize vibrations, move MagAO-X to the elevator, or place it in the Aux until ready to move to the elevator.



Slowly and carefully move MagAO-X to the telescop or Aux.

Moving the Legs to the Summit

- Return to the cleanroom with the flatbed.
- Placed on 2 dollies as in the below image.



Legs on 2 dollies placed in the middle of the table under each lower long tie bar (away from basket).

- Move the legs to the flatbed and strap them down.



Legs strapped to the truck.

- Slowly drive the truck to the summit.
- Ensure that the lift gate at the summit has been adjusted for slow smooth operation downwards (slow air release) as is done for the ASM
- Back the flatbed truck up to the lift gate.
- Next carefully roll legs onto the lift gate.
- Remove dollies from legs before going into elevator

Install MagAO-X On The Platform

- Position the telescope so that the elevator has access to the NASE platform
- Raise MagAO-X and the legs to the NASE platform
- Move the legs onto the platform, as close as possible to the telescope leaving room for lifting the instrument into place.
- Align the legs in X (parallel to the guider ring) so that only a push in Y (towards the telescope) is needed (best effort)
- Set the alignment pin system on the legs for receiving the instrument, and ensure that the platters are centered on the legs.
- Attach the load spreader using the wire harnesses as above
- Install the triangle stabilizing ropes
- Position a person at each corner of the instrument to stabilize it

- Lift the instrument with cart [1920 lbs] into position over the legs
- While keeping the instrument level, very slowly lower it into position using the alignment pins.
- If one side touches first platters will move and repeat last few steps until platters are centered and pins are centered (UPDATE once/if Jack screws arrive)
- Once on the legs, unload the crane to the cart weight of **320 lbs** and re-position the load spreader for the cart



Adjusting the load spreader for the cart.

- Now pre-load the crane again to 320 lbs.
- With the crane supporting the cart weight, remove the 4 bolts attaching the cart to the table
- Lower the cart to the floor, and detach from the crane.
- Stow the crane and handling gear
- Remove the long side of the cart on the telescope side (4 bolts), and wheel the remaining pieces out from under table. Reassemble the cart and remove to the Aux. Bldg.
- Conduct the daytime alignment procedure per alignment plan

Transport Electronics

- remove the earthquake bar
- Use the lift gate to move the electronics rack onto a truck (either the flatbed or a pickup)

- place foam between the rack side and the truck to protect cable connectors



The electronics rack has many delicate connectors on the side.

- strap the rack securely to the truck



The rack on a truck for transport.

- drive the truck to the summit
- unload the rack using the lift gate, and move to the platform on the elevator

Install Electronics and Cable

- Place the rack in position next to the instrument. Leave enough room behind it so that the door can open.
- Install the earthquake roll bar
- lock the wheels
- Check that the rack is stable and will not roll
- Install all cables but the 2K DM cables

Install AOC in Control Room

- Move the AOC, monitors and stand to the Clay control room
- Connect AOC to the telescope 200 network, and to the “VisAO” port for the internal 192 network.
- Connect the “VisAO” cable behind the MagAO rack in the equipment room directly to the media converter for VisAO.
- On the platform connect the instrument lan to the VisAO fiber with a media converter

- Power on the AOC
- Conduct function checks of everything but the 2K DM

Cable the 2K

- Cable the 2K DM following procedure

See MagAO-X PSR Document [📄 2.2 Deformable Mirrors](#)

- Check actuator functionality, following [this notebook on the RTC](#)